



Controls available

NADAL® C. difficile Tests

Rapid tests for the qualitative detection of *Clostridioides difficile* toxin A and toxin B or *Clostridioides difficile* glutamate dehydrogenase in human faecal samples.

Sample material	Sensitivity	Specificity	Product code	Pack size
NADAL® C. difficile Toxin A/B test cassette				
Faeces	> 99 %	> 99 %	582008	10 tests
NADAL® C. difficile GDH Ag test cassette				
Faeces	> 99 %	> 99 %	582004	10 tests
NADAL® C. difficile Toxin A/B + GDH test cassette				
Faeces	> 99 %	> 99 %	582022N-10	10 tests

Product features:

- ✓ Easy to use
- ✓ High sensitivity and specificity
- ✓ Fast and reliable test results in just 10 minutes
- ✓ Can be carried out using solid and fluid faecal samples
- ✓ Store test refrigerated or at room temperature (2-30 °C)



Kit contains:

- Test cassette
- Sample collection tube with buffer
- Package insert

Our service team would be happy to answer any questions you may have by telephone on: +49 941 29010-0 or +49 2841 99 820-0 and by email at: vertrieb@nal-vonminden.de

NADAL® C. difficile Test



Clostridioides difficile infection

The bacterium *Clostridioides difficile* (*C. difficile*) is usually transmitted via the faecal-oral route. Its two toxins, toxin A and toxin B, are responsible for its symptoms, which can range from mild, watery diarrhoea to pseudomembranous colitis or even toxic megacolon and can potentially lead to death. *Clostridioides difficile* glutamate dehydrogenase (*C. difficile* GDH) is an enzyme which is produced in large quantities by all toxigenic and non-toxigenic strains, making it an excellent biomarker.


NADAL® C. difficile rapid tests

The NADAL® C. difficile rapid tests are immunochromatographic assays in lateral flow format for the qualitative detection of *C. difficile* toxin A and toxin B or *C. difficile* GDH in human faecal samples.


The test is best suited for a diagnosis of *C. difficile* infections in hospitals and laboratories. With a short test time of just 10 minutes, as well as clear, easy-to-read results, these tests are the ideal aid in the diagnosis of a *C. difficile* infection.

We help you care!

Test procedure*

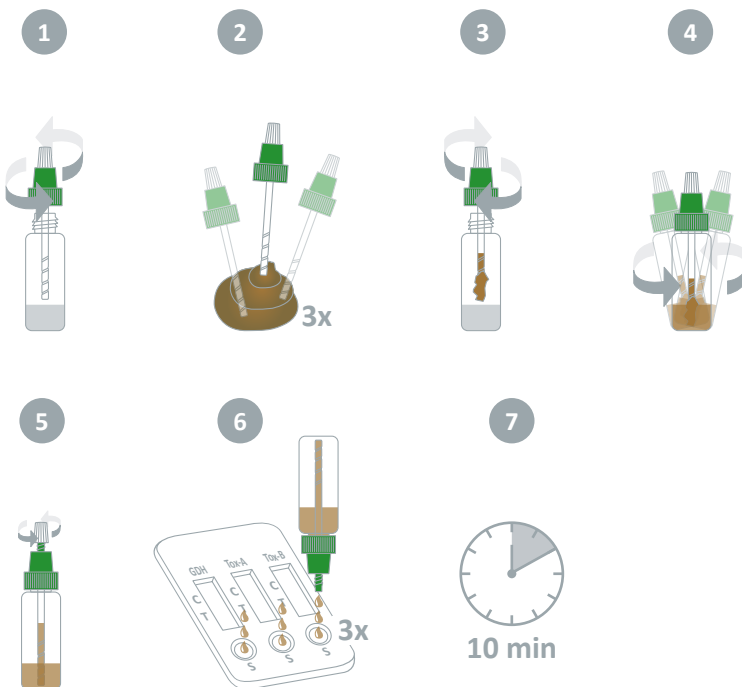

Result interpretation

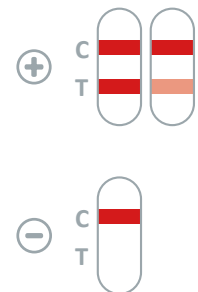

A detailed description can be found in the package insert.

*There is a slight difference in the test procedures between the NADAL® C. difficile GDH Ag Test and NADAL® C. difficile Toxin A/B Test.

